

# WATCHHALER™ - TURNS INHALATION INTO CHILD'S PLAY



**BETTER DESIGN**  
FOR BETTER COMPLIANCE

**VISIBLE FEEDBACK**  
TO PARENTS

**FLOW AND**  
VOLUME CONTROL

Current spacers are limited by a sterile design approach, while children live in a different world. Activaero's Watchhaler wants to be a "helping friend" for asthmatic children.

**MARKETING AND**  
DISTRIBUTION PARTNERS  
WELCOME!

- For further information on
- technology
  - partnerships
  - cooperations

please visit us at  
**STAND A37**



**WATCHHALER™**



**ACTIVAERO**  
TECHNOLOGIES

# WATCHHALER™ - TURNS INHALATION INTO CHILD'S PLAY

## BETTER DESIGN AND BETTER TECHNOLOGY




The Watchhaler can easily be disassembled and enables a thorough and quick cleaning process.

The new Watchhaler contains an inner balloon that controls the inhalation volume. Inspiratory flow is limited to 150 ml/s by Activaero's patented flow limitation valve.

## SYSTEM COMPARISON

	MDI	SPACER	WATCHHALER
Application of an efficient Dose	✓	✓	✓
Safety	✓	✓	✓
Easy coordination and handling		✓	✓
Safe, reliable and constant Dose			✓
Parents can easily control the inhalation procedure			✓
electro-static charge	N/A		✓
Child-adequate Design			✓

## EASY HANDLING VISIBLE CONTROL

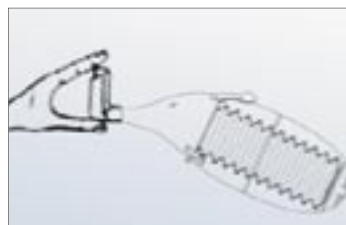

A single puff of an MDI is released into the Watchhaler.



The child and its attendant see the deflation of the balloon during inhalation.



The fully deflated balloon signals a successful inhalation.

## OUR OFFER

Watchhaler is protected by several international patents. Activaero is looking for marketing and distribution partners with a focus on asthma treatment in children.

Interested in Watchhaler?  
More information are available at:

**Activaero GmbH**  
Wohraer Str. 37  
35285 Gemünden  
Germany

email: [watchhaler@activaero.de](mailto:watchhaler@activaero.de)  
internet: [www.activaero.de](http://www.activaero.de)



**WATCHHALER™**



**ACTIVAERO**  
TECHNOLOGIES